


Plantweb Optics and AMS Machine Works Software Instructions

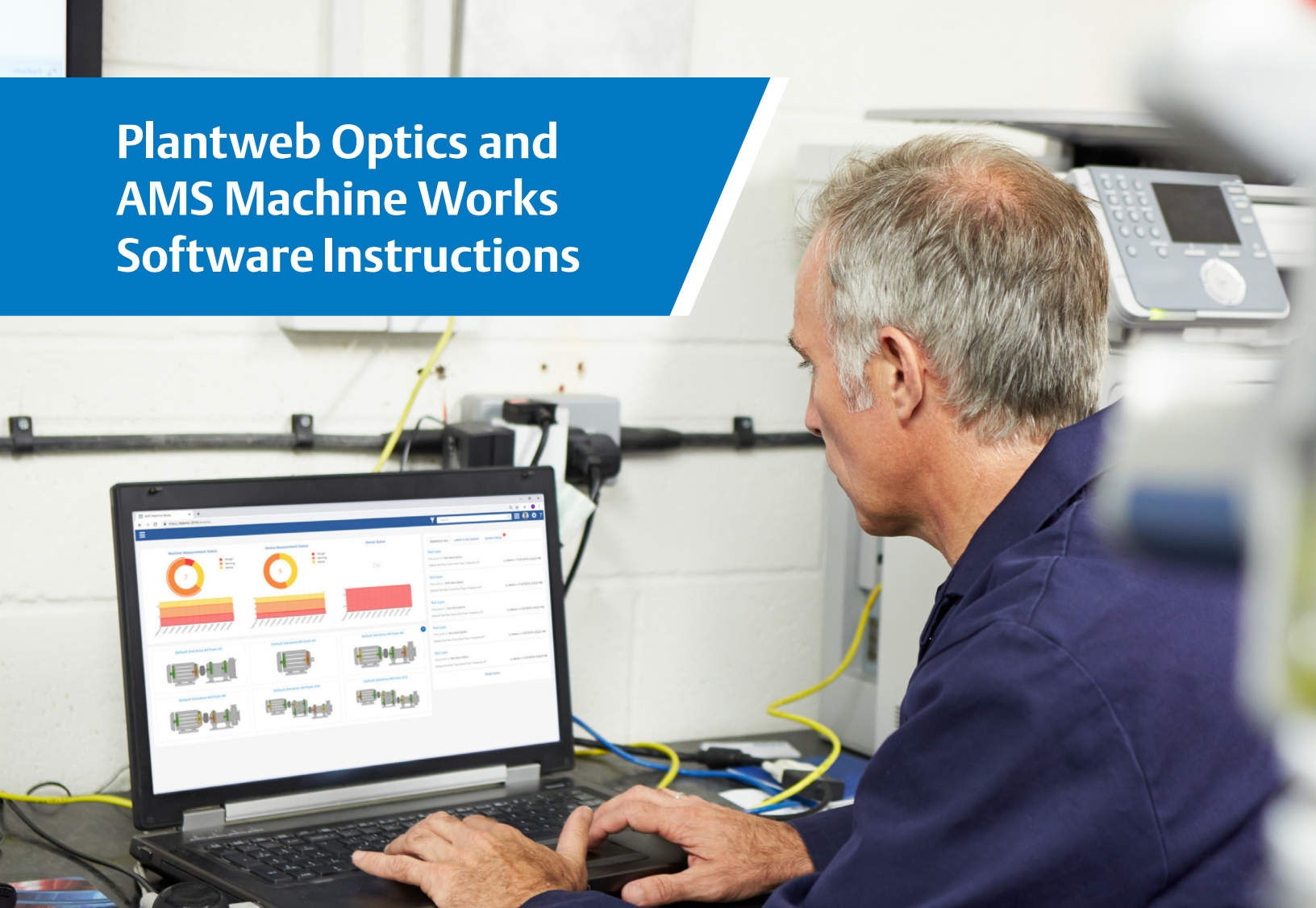

To register your Plantweb Optics or AMS Machine Works software:

The machine fingerprint, or lock code, of the server where Plantweb Optics is installed must be retrieved before you can generate a license file. The echoid tool found in the Plantweb Optics installation zip file is used to retrieve the machine fingerprint or lock code.

Prerequisites

- A48OPTICS-SYSTEM0.Plantweb_Optics.1.6.X.X.zip file.

Procedure

1. Extract the A48OPTICS-SYSTEM0.Plantweb_Optics.1.6.X.X.zip file to the Plantweb Optics server.
Note: Extract the zip file on a root directory.
For example, drive C.
2. Start a command prompt with administrator privileges.
3. Change directories to the following directory in the Plantweb Optics installer zip file: A48OPTICS-SYSTEM0.Plantweb_Optics.1.6.X.X\Setup_Support\Gemalto
4. Run echoid.exe -d -p from the command prompt window to generate your locking code.